

How to treat your LARGEST patients as well as seniors with the new Shuttle Recovery



The Shuttle Recovery has been specifically designed to accommodate your largest bariatric and physical therapy patients while providing the comfort and confidence your smaller senior clients require.

Extra Wide, Extra Strong Carriage

Start with an oversized carriage that is a full two feet wide – our largest ever. The bariatric patient can be completely confident there is ample room. Because of its extra width, seniors also feel very safe in lying down for treatment. The carriage has a hard foam layer topped with soft foam then covered in a heavy duty vinyl so patients are also relaxed and comfortable during their treatment.

The carriage is supported by six polyurethane wheels which allow the Recovery to accommodate a patient weight up to 500 pounds. These wheels run through aluminum tracks constructed with a special box configuration – unique to Shuttle products – that provides additional strength while ensuring a smooth run. The Recovery gives you strength as well as a smooth and very quiet operation.

Adjusts Around the Patient

Once on, you can adjust the Recovery to fit the patient. ROM control allows you to adjust the relative position of the carriage to the kickplate so you can control the patient flexion angles. Elasticord resistance is all engaged at the head-end so the patient does not have to get on and off to change their resistance.

The elasticords allow you to provide progressive treatments from 12½ lbs. all the way up to 300 lbs. at full extension. Use it for non-weightbearing protocols, leg press exercises and plyometrics. Yes, bariatrics and seniors can do plyometrics on the Recovery (and they couldn't get this kind of a workout in any other environment without significant impact to their joints and spine).

Affordable, Smaller Footprint

So, the Recovery can certainly accommodate your large patients – and it can suit all your other patients at the same time. What's more, we have designed it to function within the increasing space and fiscal demands of the physical therapy environment.

The Recovery is our smallest “big” machine measuring just eight feet in length. Even with its added strength, the Recovery weighs in at less than 200 pounds. Finally, and perhaps best of all, the Recovery starts at an affordable price.

To learn more about how the versatile Shuttle Recovery can help you welcome all shapes and sizes of patients, please call **1-800-334-5633** or visit our website at **ShuttleSystems.com**.

Join us online:



Contemporary Design Company
SHUTTLE RECOVERY

SHUTTLE® RECOVERY

“Treat Across the Age and Ability Spectrum”

Padded Kickplate

Oversized for a variety of activities, the 24” x 30” kickplate has a removable rubber pad with soft foam for joint protection.



**Included with Recovery SR model or can be purchased separately.*

Recovery SR Foam Wedge Backrest*


Covered with heavy-duty, soft vinyl for patient comfort and confidence. Is easy to clean and ultra soft to the touch. Filled with high compression urethane foam that will keep its shape over repeated use.

Oversized Carriage

Measuring 2’W x4’L at a height of 21”,with a capacity of 500 pounds, it accommodates patients of all sizes with comfort and confidence.

Head-end Resistance System

Lanyards at the head-end makes it easy for practitioners to change resistance loads during treatment without requiring the patient to get off and on the equipment.

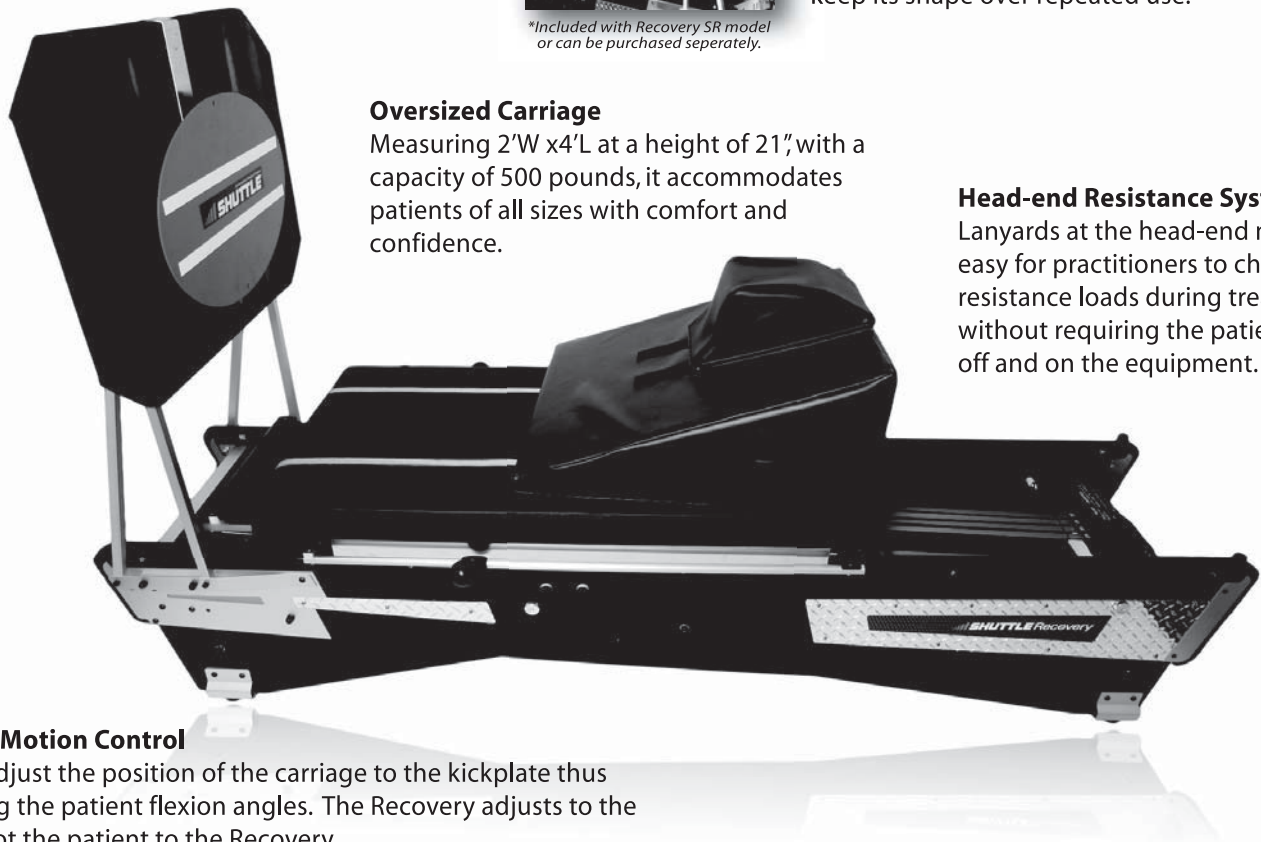

Range of Motion Control

You can adjust the position of the carriage to the kickplate thus controlling the patient flexion angles. The Recovery adjusts to the patient, not the patient to the Recovery.

Features:	Recovery Standard 5200	Recovery Sport 5300	Recovery Performance 5400	Recovery Senior 5500
Progressive Elastic Resistance	12.5 - 150 lbs. @ full extension	12.5 - 300 lbs. @ full extension	12.5 - 300 lbs. @ full extension	12.5 - 300 lbs. @ full extension
Widebody Carriage (500 lb. capacity)	24” x 48” 6 wheel	24” x 48” 6 wheel	24” x 48” 6 wheel	24” x 48” 6 wheel
Headrest	Low Profile Headrest*	Universal Headrest	Universal Headrest	Universal Headrest
Plyometric Steel Brace		Steel Cross Brace	Steel Cross Brace	Steel Cross Brace
PNF Pulley System with Towers			x	
24” Horizontal Structure				x
Accessories				
Shuttle Wobble Board				x
Foam Wedge Backrest				x
Kickplate Foot Supports				x
Graphic Kickplate Covers Available				
Treatment Table Extension (24” x 24”)				

* Not shown in image

ShuttleRecovery Specifications

Overall Footprint: 31” x 98” * Backrest Dimensions 24” x 48” * Resistance Load: 12.5 - 300lbs.
 Carriage Capacity: 500 lbs. with 6 wheels * Warranties: Structural - 2 years * Upholstery - 1 year * Parts 1 year
 1-800-334-5633 * www.ShuttleSystems.com